



PRODUCT CATALOG



MADE IN THE USA 



CHECK OUT
OUR NEW
SIZES!

MORE PROTECTION. LESS MOUTHGUARD.™

sisuguard.com



MOUTH GUARDS

UP TO \$50,000 DENTAL WARRANTY ON ALL OF OUR MOUTHGUARDS
CUSTOM FIT AND REMOLDABLE UP TO 20 TIMES

AERO FOR TEAM SPORTS

At a lean 1.6mm, the ultralight design of the Aero Mouthguard is 50% thinner than conventional sports mouth guards, providing maximum comfort and convenience. Combined with its quick, custom and remoldable fit, you'll entirely forget it's there.



Classic Shape



1.6mm thin

Made in the USA

SMALL

Up to 5' (150 cm)
For smaller mouths and children

Available in 13 colors:



MEDIUM

5' - 6' (150 cm - 180 cm)
For average mouths

Available in 16 colors:



LARGE

Over 6' (180 cm)
For larger mouths

Available in 6 colors:



MAX FOR HIGH-IMPACT SPORTS

At 2.4mm thin, the SISU Max Mouthguard offers the most protection of any SISU guard, but is still 30% thinner than conventional sports mouth guards. It's custom-fitted so you'll be able to talk, breathe and drink better than ever before. It's really the best of both worlds.

Available in 8 colors:



Classic Shape



2.4mm thin

ONE SIZE Over 5' (150 cm)

For average and larger mouths

Made in the USA



GO FOR ESSENTIAL PROTECTION

Smaller in size and price, the all new SISU GO Mouthguard was made to be both discreet and comfortable while still providing the necessary canine to canine coverage to protect where most dental injuries happen. The SISU protection athletes expect, at a more affordable price point.

Available in 5 colors:



Low-Profile Shape



1.6mm thin


ONE SIZE Over 4' (120 cm)

One size fits most

Made in the USA

SENSE ELECTRONIC HEAD IMPACT TRACKER

With the SISU Sense™ Mouthguard, we've taken our mission to protect the athlete to the next level. By infusing our HitSense™ Technology into our industry-leading mouthguard, SISU Sense measures every head impact an athlete absorbs, then transmits the data to our proprietary app. **The result: protection on the field today, and tomorrow, that makes sense.**

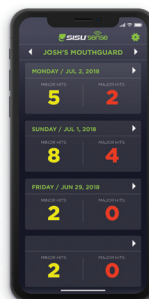

Classic Shape



2.4mm thin

Available in 3 colors:



ONE SIZE Over 5' (150 cm)

For average and larger mouths

Made in the USA

CASE

Lightweight and durable
Pin hinge
Extra storage



MOUTHGUARD STORAGE

GUARD SOLD SEPARATELY

SISU FRESH

Cleans and freshens
Cinnamint flavor
Convenient, ready to use spray bottle
Click-on cap to prevent leaking
About 450 sprays per bottle
2 oz



CLEANS + FRESHENS

spray. compete. repeat.

Made in the USA

SISU HEAT PACK

Convenient, portable heating solution
Ready to mold in 3 minutes
No water or electricity needed
Works with flat mouthguards only
1-time use




PORTABLE HEATING

Mold on-the-go!

Made in the USA

FLAT2FIT™ TECHNOLOGY

SISU Mouthguards use Flat2Fit technology to create a completely custom fit to your mouth. The material comes flat in the packaging, and within minutes, transforms into a perfectly fit mouth guard ready to protect you from whatever comes your way.

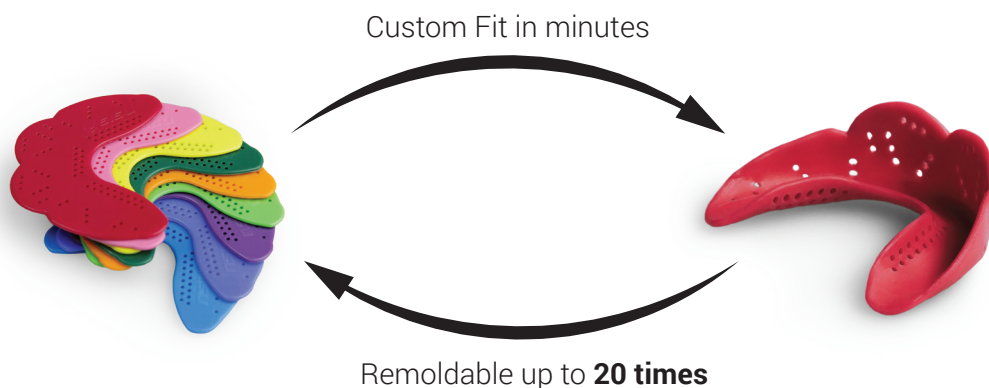




POWERED BY DIFFUSIXTM TECHNOLOGY

THE SISU ADVANTAGE:


SISU Mouthguards are made from a non-compressible, perforated material that distributes impact force over the entire surface of each mouthguard. With their custom, remoldable fit and unmatched strength, SISU Mouthguards are able to better protect your teeth while providing you the freedom to easily talk, breathe and drink.

Conventional Mouthguard
Impact Force Dispersion

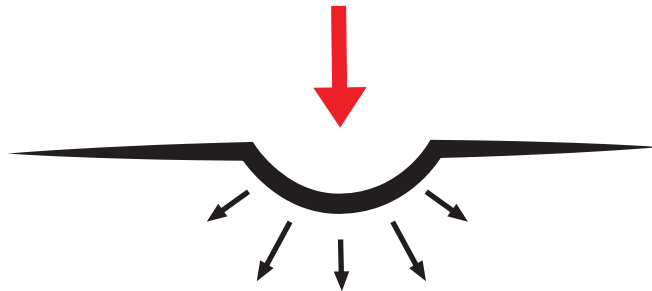

Conventional mouthguards are made from a material that is too soft to fully protect your teeth from harm.

SISU
SISU Mouthguards
Impact Force Dispersion

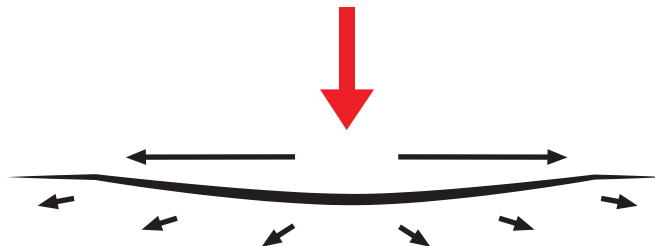

With our one-of-a-kind DiffusixTM Technology, SISU Mouthguards are made from a material that absorbs and diffuses impact throughout the entire mouthguard, meaning that much less impact force is transferred to your teeth.






QUICK-FIT GUIDE

Notice:

If you have braces, **DO NOT** attempt to fit your SISU Mouthguard at home. Instead, contact your orthodontist (see below for more info).

If you don't have braces, let's get started.

- 1** Wash your hands and the SISU Mouthguard with soap and water. Heat water in a small bowl until hot to the touch (160°F / 71°C). Place your SISU Mouthguard into the hot water until soft and pliable, approx. 30-60 seconds (use a non-plastic utensil to determine softness).
- 2** Once soft, use the utensil to carefully lift the SISU Mouthguard out of the water. **DO NOT let the softened mouthguard fold on itself – the material will glue together, making it difficult to reshape.**
- 3** Balancing the pliable guard with your thumbs, lean in toward the mirror and open your mouth just wide enough to insert the guard. Insert the SISU Mouthguard and close down gently on the bite pad in the middle of the guard.
- 4** Hold the mouthguard in place using your teeth and tongue. With your fingers, fold the guard straight up. Be sure to mold the guard up to the gum line on the front and sides of your teeth to create a solid shield. **DO NOT** stretch the guard.
- 5** Once the SISU Mouthguard is in the proper position, close your lips over the guard and gently create suction until the material has hardened to a snug, perfect fit. Leave the guard in your mouth for 3-5 minutes to allow time to set.

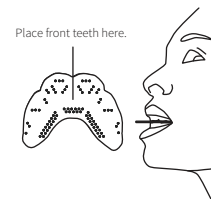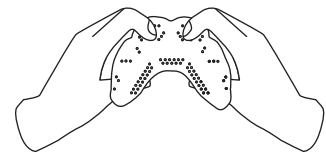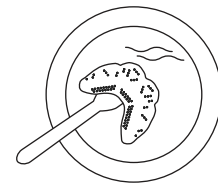

GOT BRACES?

If you have braces or are undergoing other orthodontic treatment, **DO NOT** attempt to fit your SISU Mouthguard at home. Instead have your dentist or orthodontist fit your SISU Mouthguard for you! He/she will use the appropriate barriers between the guard and braces so that the material, while pliable, doesn't bond to your brackets.

We also recommend that no matter your child's age, all athletes use an Aero Medium or Aero Large SISU Mouthguard whenever wearing braces in order to ensure a proper fit.

Visit sisuguard.com/fitting-with-braces for more information.



THE OFFICIAL MOUTHGUARD OF

MOUTHGUARD COMPARISON



AERO

MAX

GO

SENSE

Sizes	Small	Medium	Large	One Size	One Size	One Size
Athlete Height	Up to 5' (150 cm)	5' - 6' (150 cm - 180 cm)	Over 6' (180 cm)	Over 5' (150 cm)	Over 4' (120 cm)	Over 5' (150 cm)
Works with braces	No (Use Medium or Large)	Yes	Yes	Yes	Yes	Yes
Thickness	1.6 mm			2.4 mm	1.6 mm	2.4 mm
Dental warranty	Up to \$50,000			Up to \$50,000	Up to \$50,000	Up to \$50,000
Uses	All sports			High-impact sports	Essential Protection	Head Impact Tracking
Best for	Lacrosse Field Hockey Ice Hockey Football Basketball Roller Derby Team Sports And more!			Rugby MMA Boxing Karate/BJJ Hockey Football Impact Sports And more!	Basketball Soccer Cycling Water Polo Baseball Volleyball Action Sports And more!	Rugby MMA Lacrosse Roller Derby Hockey Football Impact Sports And more!
MSRP	\$24.99			\$34.99	\$11.99	\$99.99

Research suggests the taller you are, the larger your dental arch. And thus, a bigger mouthguard is required as height increases. Our new and improved sizing system uses the findings of extensive research and testing to suggest the proper guard for you based on how you measure up.

*SISU Mouthguards are BPA, latex, PVC, and phthalate free, and FDA compliant.

CONTACT INFORMATION

Questions? Please contact our team at 800.444.0570
or send an email to customerservice@sisuguard.com

sisuguard.com